



forza Heavy-Load Bed
www.muehle-mueller.de

forza heavy-load bed

Nursing and hospital version

Nursing and care of overweight or adipose patients often is connected to increased treatment risk and high effort for the nursing staff.

Overweight patients – with a weight of 150 kg and more – are often not the exception in nursing work. This brings special challenges for nursing staff as well, since standard nursing beds or aids often do not have enough load capacity.

3 Model Versions

The forza nursing bed series – with the models forza basic, forza trend and forza aktiv – facilitates everyday work for nursing staff and family members acting as nurses.

Extremely stable construction – load up to 500 kg

With a load capacity of up to 500 kg and the many adjustment options of a “normal” nursing bed, forza permits perfect positioning and freedom of movement for heavy persons. With its stable construction, it is ready for any demands of everyday nursing work.

forza in 3 different versions

- Standard model **forza basic** with a height adjustment of 40 - 80 cm
- **forza trend** with head and foot lowering
- Top model **forza aktiv** as sitting and rising bed with exit on the foot end

Comfort

- High lying comfort with 4-fold separated steel grid lying area and special mattresses for overweight /adipose persons requiring nursing

Side grid

- Continuous side grids
- Optionally with 2-part side grids

Function

- Available in versions with a load capacity up to 250 kg, 350 kg and 500 kg
- Height adjustment 40 - 80 cm
- Trendelenburg / Anti-Trendelenburg
- Getting out at the foot end
- Mains release
- Single or 4-wheel central parking
- Hand switch can be locked
- Complete bottom clearance



Areas of use

- In nursing at home, in seniors' residences and nursing homes, hospices, rehabilitation facilities, hospitals and retirement homes



Design

- Comfortable and beautifully shaped head and foot parts

Individuality

- Other wood patterns on request
- Metal coating in other RAL colours
- Lying surface widths 115 cm or 120 cm (140 cm only forza basic and forza trend)
- Bed can be extended to 220 cm (only forza basic)

forza aktiv



forza basic

Electrical adjustment

- Of the height
- Of the back, thigh and lower leg rests



forza trend

Electrical adjustment

- Of the height
- Of the back, thigh and lower leg rest
- Trendelenburg / Antitrendelenburg



forza aktiv

Electrical adjustment

- Of the height
- Of the back, thigh and lower leg rests Trendelenburg / Antitrendelenburg
- Sitting and rising positions

forza basic heavy-load bed

The strong basic model

forza basic is the standard model of our forza nursing bed series. The heavy-load bed has an electrically adjustable back, thigh and lower leg rest and can be adjusted from 40 - 80 cm in height.

- Electrically adjustable back, thigh and lower leg rest via hand switch
- Height adjustment of 40 to 80 cm
- 4-fold separated steel grid lying surface
- Ground clearance for patient lifter
- Available with continuous or 2-part side grids
- Carriers up to 250 kg, 350 kg or even 500 kg
- Stable and vibration-free
- Standard pattern: light beech other patterns on request



Setting the lying surface low



Raising the lying surface



Head and foot section adjustment




forza trend heavy-load bed

Lowering of foot or head end (Trendelenburg / Antitrendelenburg)

The forza trend has a Trendelenburg function (lowering head/foot end) and convinces with its stable construction. It has all the standard functions of a nursing bed and can be individually adjusted in height between 40 and 80 cm. Ideal for everyday nursing of adipose and overweight persons.

- Electrically adjustable back, thigh and lower leg rest via hand switch
- Lowering head/foot end (Trendelenburg/Anti-Trendelenburg lockable)
- Height adjustment from 40 to 80 cm
- Available with continuous or 2-part side grids



Lowering the lying surface



Raising the lying surface



Lowering foot end (Trendelenburg)



Lowering head end (Anti-Trendelenburg)



forza aktiv heavy-load bed

The sitting and rising bed with its exit at the foot end

With its exit on the foot end, the heavy-load bed forza simplifies getting into and out of the nursing bed. Many other adjustment options support the nursing staff and family performing nursing tasks especially with the different types of patient nursing, particularly facilitating care in the genital area, of the feet and the entire body. The options for perfect positioning and movement of the patient permit comfortable sitting and rising positions.

- Electrically adjustable back, thigh and lower leg rest via hand switch
- Height adjustment from 40 to 80 cm
- Sitting and rising function electrically adjustable via hand switch
- Beautifully shaped handholds for getting into and out of the heavy-load bed easily
- Removable, beautifully shaped foot bed
- 4-fold steel grid lying surface
- Available with continuous or 2-part side grids
- Carriers up to 250 kg, 350 kg or even 500 kg
- Stable and vibration-free

Functions



Lowering the lying surface/lying position



Raising the lying surface



Resting position



Sitting and nursing position



Sitting and exit position



forza aktiv mit durchgehenden Seitengittern



forza klinik

All benefits of the forza hospital version

The forza heavy-load bed is also available in a hospital version.


All properties and options of the basic function apply here as well.

In hard everyday use in the hospital, stability, load capacity and function are needed in any situation.

The heavy-load bed forza by Mühle-Müller easily meets these requirements and turns out to be a reliable companion for patient and nursing staff.



forza Wood Patterns

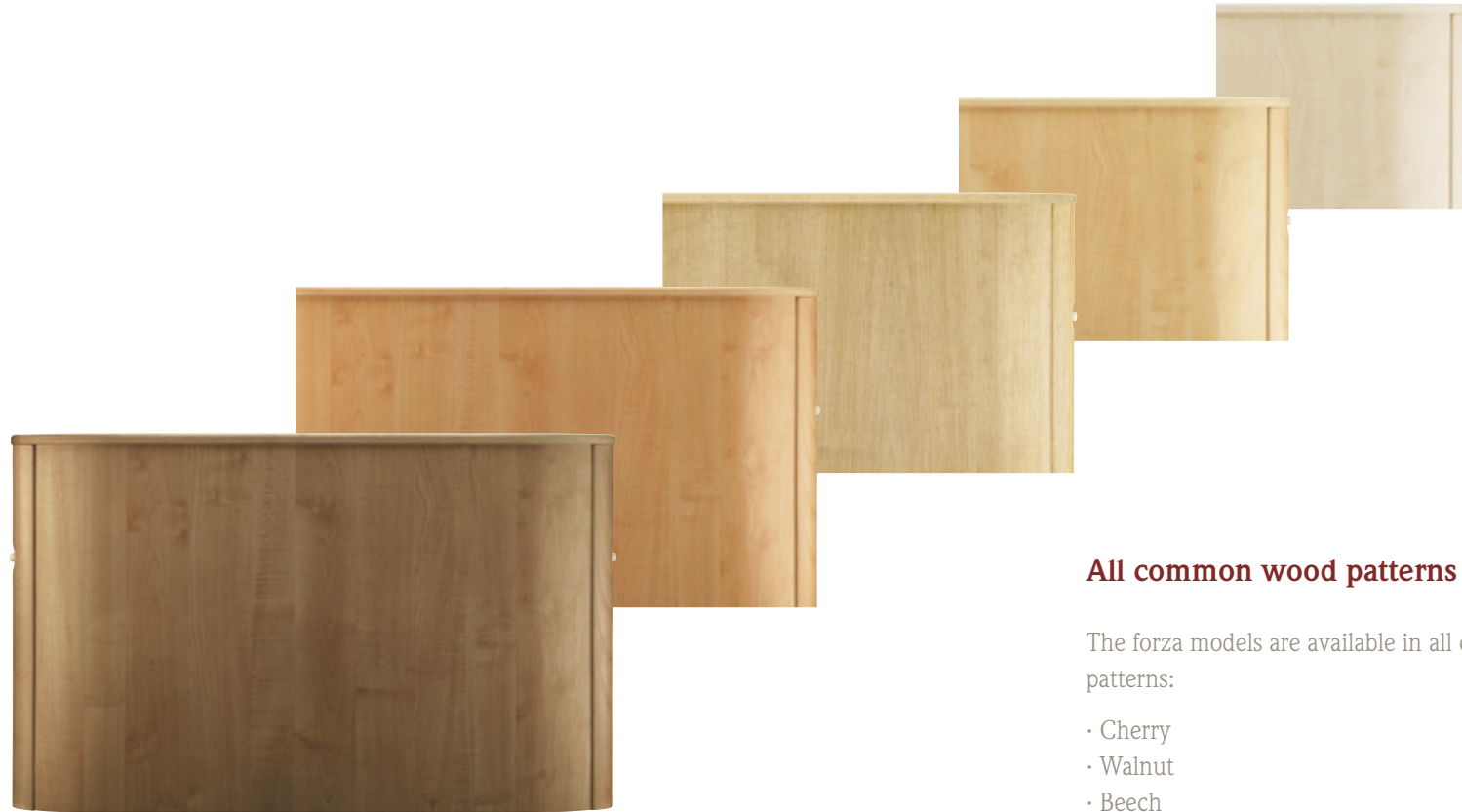

All common wood patterns

The forza models are available in all common wood patterns:

- Cherry
- Walnut
- Beech
- Oak
- Birch
- And many others



Production

Werner Mühle Metallverarbeitung
Boschstrasse 6
35799 Merenberg

Phone: +49 (0) 6471 - 5065-0
Telefax: +49 (0) 6471 - 5065-13

E-Mail: info@muehle-mueller.de
Internet: www.muehle-mueller.de



Sales

m2 Handels- u. Vertriebs GmbH
Oststrasse 26
51766 Engelskirchen - Ränderoth

Phone: +49 (0) 2263 - 9617-0

E-Mail: info@muehle-mueller.de
Internet: www.muehle-mueller.de